

2015


## Workflow Modifications and Addition of MALDI-TOF Technology Significantly Improved Turn-Around-Time to Identification of Common Urine Pathogens

Dimath Alyemni, MD  
*Thomas Jefferson University*

Christine M. Lam  
*Thomas Jefferson University*

Amity L. Roberts, PhD, D(ABMM)  
*Thomas Jefferson University*

Follow this and additional works at: <https://jdc.jefferson.edu/pacbresidentposters>

 Part of the [Medical Anatomy Commons](#), [Medical Cell Biology Commons](#), and the [Medical Pathology Commons](#)

[Let us know how access to this document benefits you](#)

### Recommended Citation

Alyemni, MD, Dimath; Lam, Christine M.; and Roberts, PhD, D(ABMM), Amity L., "Workflow Modifications and Addition of MALDI-TOF Technology Significantly Improved Turn-Around-Time to Identification of Common Urine Pathogens" (2015). *Department of Pathology, Anatomy, and Cell Biology Resident's Posters*. Paper 14.

<https://jdc.jefferson.edu/pacbresidentposters/14>

This Article is brought to you for free and open access by the Jefferson Digital Commons. The Jefferson Digital Commons is a service of Thomas Jefferson University's [Center for Teaching and Learning \(CTL\)](#). The Commons is a showcase for Jefferson books and journals, peer-reviewed scholarly publications, unique historical collections from the University archives, and teaching tools. The Jefferson Digital Commons allows researchers and interested readers anywhere in the world to learn about and keep up to date with Jefferson scholarship. This article has been accepted for inclusion in Department of Pathology, Anatomy, and Cell Biology Resident's Posters by an authorized administrator of the Jefferson Digital Commons. For more information, please contact: [JeffersonDigitalCommons@jefferson.edu](mailto:JeffersonDigitalCommons@jefferson.edu).



Abstract

**Background:**  
In order to improve the identification of common aerobic urine cultures as well as antimicrobial susceptibility testing (AST) setup at an Academic Medical Center, work-flow modifications and MALDI-TOF technology were incorporated. Previously, the majority of species identification was achieved with conventional identification/antimicrobial susceptibility combo panels. All urine cultures, regardless of laboratory receipt time, were previously read once per day on 1<sup>st</sup> shift.

**Methods:**  
The initial workflow modification involved addition of a 2nd shift urine culture reading. Urine specimens received from 8:00 AM to 4:00 PM were read on 1<sup>st</sup> shift, while urine specimens received from 4:00 PM to 8:00 AM were read on 2nd shift.


Additionally, urine cultures were sorted into categories: no growth (NG) at 24 hours, no growth at <24 hours, single colonies of growth, multiple colonies of growth, and potential contaminants. No growth cultures were signed out at 24 hours. No growth cultures at < 24 hours were reincubated to be read on subsequent shift. Cultures with growth were set aside as either single colony types or multiple colony types. Cultures of probable contaminants were signed out.

Once cultures were sorted, the isolated colonies underwent MALDI-TOF analysis (Bruker) and antimicrobial susceptibility testing (AST) as appropriate. Individual technologists setup the MALDI-TOF target plate map and spotted the associated target plate. AST was setup at the same time. The MALDI-TOF was then operated by a central technologist and results reported by the original technologist reading the culture.

**Results:**  
Retrospective pre-workflow (September-November 2013) and post-workflow (May, June, October 2014) modification turn-around-times were compared for 16 commonly isolated pathogens. These pathogens consisted of common urine pathogens as noted in Table 1. *Staphylococcus aureus* was previously identified in our laboratory by a positive coagulase test and not included in this analysis. The average turn-around-times, standard deviations and the p-values for each organism are indicated in Table 1.

**Conclusion:**  
Converting from conventional identification methods to MALDI-TOF, in conjunction with workflow modifications such as a 2<sup>nd</sup> culture reading, notably improved urine culture turn-around-time for identification and AST.

Laboratory Process

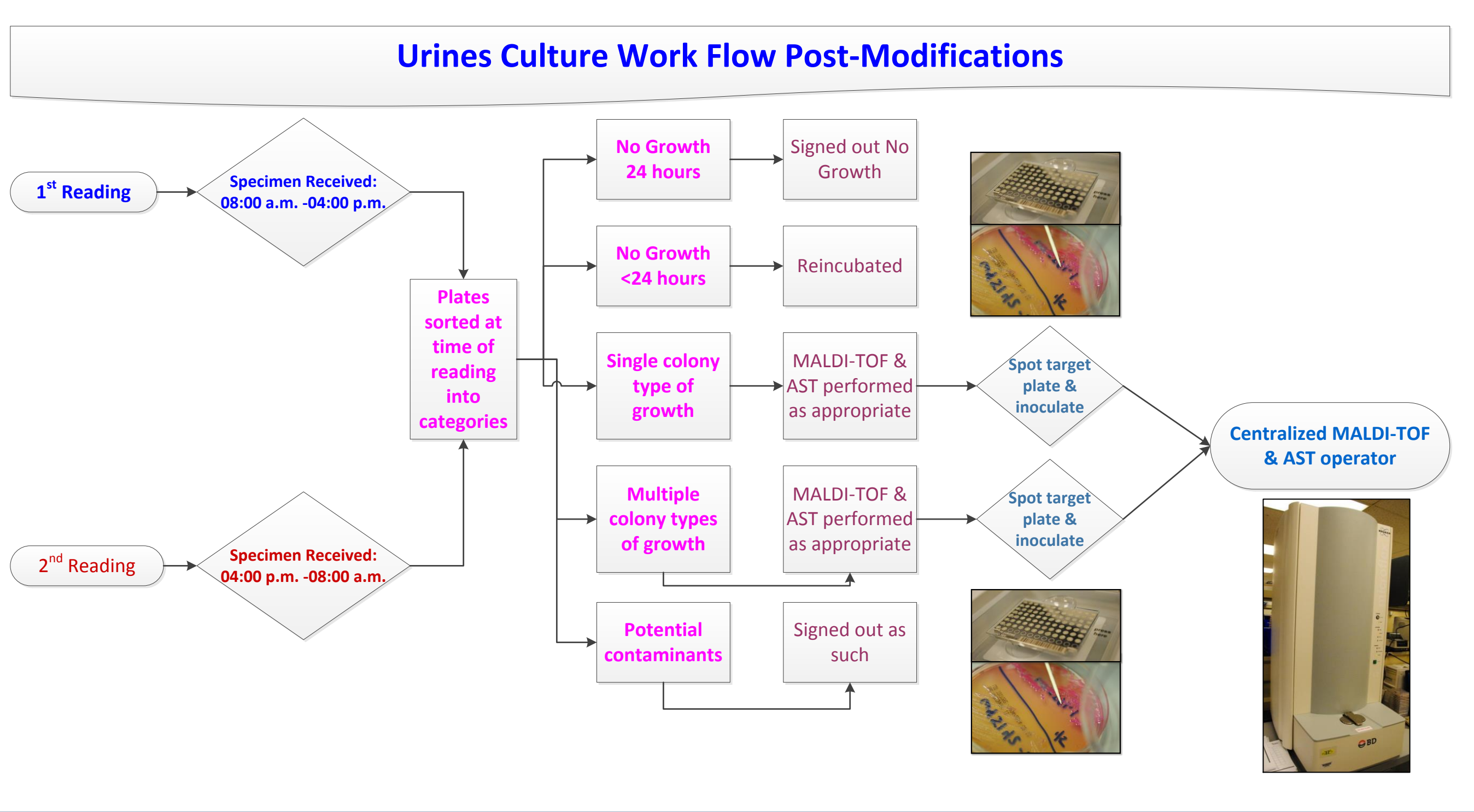

Motivation for implementation:

- Surgical and Oncology Clinicians were dissatisfied with the urine culture turn-around-time, especially in cases of 24 hour negative cultures as this impacted discharge for high risk patients
  - The solution from the lab perspective was to incorporate an afternoon reading for specimens received later on the prior day, allowing 24 hour negative cultures to be signed out before accrual of an additional inpatient day
- Further improvements were sought through incorporation of MALDI-TOF technology to allow faster identification times.
  - This also allowed clinicians to apply more effective empiric antimicrobial therapy based on organism identification.

Table 1. The average turn-around-times for common urine pathogens

Organism	Pre-Workflow Modifications (Sept., Oct., Nov., 2013)		Post-Workflow Modifications (May, June, Oct. 2014)		P-value (<0.05 significant)
	Average hrs (N)	STDev (±)	Average hrs (N)	STDev (±)	
<i>Acinetobacter baumannii/ calcoaceticus complex</i>	62.42 (12)	24.36	54.65 (11)	18.56	0.398
<i>Citrobacter freundii complex</i>	73.93 (22)	22.36	39.24 (30)	48.52	5.792x10 <sup>-07</sup>
<i>Citrobacter koserii</i>	62.81 (26)	21.70	29.46 (27)	45.52	2.559x10 <sup>-08</sup>
<i>Enterobacter aerogenes</i>	62.73 (32)	17.14	27.87 (21)	8.49	2.194x10 <sup>-10</sup>
<i>Enterobacter cloacae complex</i>	62.9 (45)	23.80	40.36 (41)	22.65	2.128x10 <sup>-05</sup>
<i>Enterococcus faecalis</i>	67.12 (246)	23.28	38.85 (164)	19.63	2.01x10 <sup>-29</sup>
<i>Enterococcus faecium</i>	91.39 (19)	25.75	57.06 (67)	29.83	2.163x10 <sup>-05</sup>
<i>Escherichia coli</i>	52.83 (1041)	15.22	30.50 (929)	12.92	4.29x10 <sup>-63</sup>
<i>Klebsiella oxytoca</i>	53.31 (28)	15.39	33.59 (15)	16.56	0.001
<i>Klebsiella pneumoniae</i>	55.93 (291)	18.02	38.52 (255)	22.43	3.50x10 <sup>-21</sup>
<i>Morganella morganii</i>	66.83 (21)	19.96	36.94 (14)	14.00	1.016x10 <sup>-05</sup>
<i>Proteus mirabilis</i>	58.74 (143)	21.19	31.05 (109)	11.54	4.01x10 <sup>-30</sup>
<i>Pseudomonas aeruginosa</i>	66.82 (155)	20.47	40.54 (88)	21.28	3.42x10 <sup>-17</sup>
<i>Serratia marcescens</i>	62.37 (11)	23.85	29.95 (16)	24.15	0.001
<i>Staphylococcus saprophyticus</i>	69.50 (21)	16.40	55.72 (19)	23.31	0.040
<i>Stenotrophomonas maltophilia</i>	85.71 (10)	20.49	58.18 (8)	46.95	0.018

Future Considerations:

- Evaluate clinical history for urines collected on inpatients to determine rate of urosepsis
  - Evaluate the difference between Length of Stay (LOS) in the pre-modification verses post-modification group
  - Trend in decreased length of stay noted, but further analysis to needed to determine root cause.
- Future algorithm, in conjunction with Infection Control, will be incorporation of a positive Urinalysis with reflex to Urine Culture and a negative Urinalysis with cancellation of Urine Culture